




PRESTO
Transportation trolleys

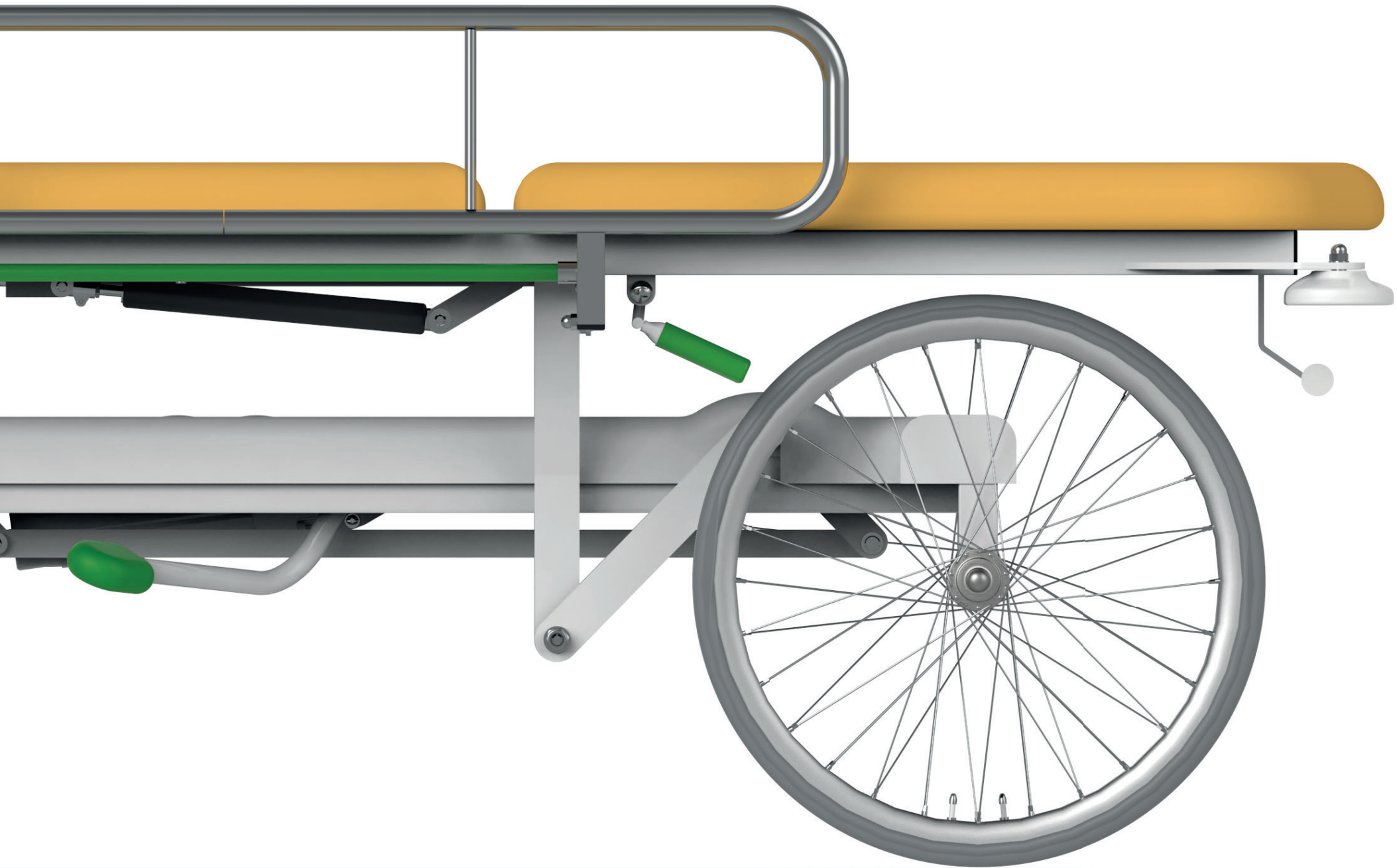

PRESTO

The PRESTO transport trolley is designed for transporting immobile patients. Ergonomic and intuitive controls allow immediate use of all functions. The large castors and lightweight design ensure high mobility and easy handling.




High load capacity

The solid and robust design with a load capacity of up to 230 kg creates a safe support for patient care.



Mattresses made of medical leatherette

The mattress platform consists of a mattress made of medical leatherette. The leatherette has an antistatic treatment that is antibacterial and has excellent antimicrobial properties.

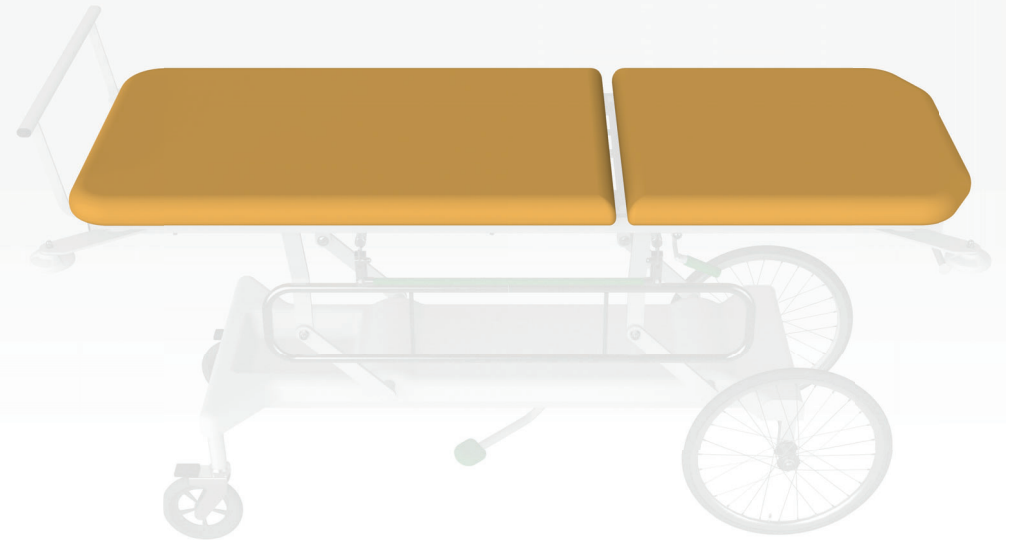



Choice of colour

The option to choose the colour design of the trolley's leatherette allows you to create an aesthetic and welcoming space that can help to improve mood, comfort and reduce stress.



Easy height adjustment

The mattress platform lift is controlled by an easy-to-use foot pedal accessible from both sides.

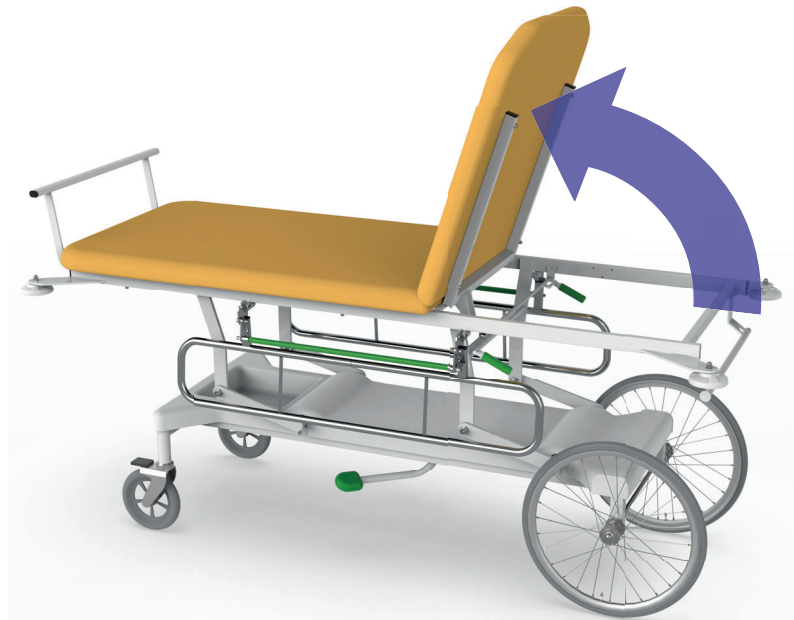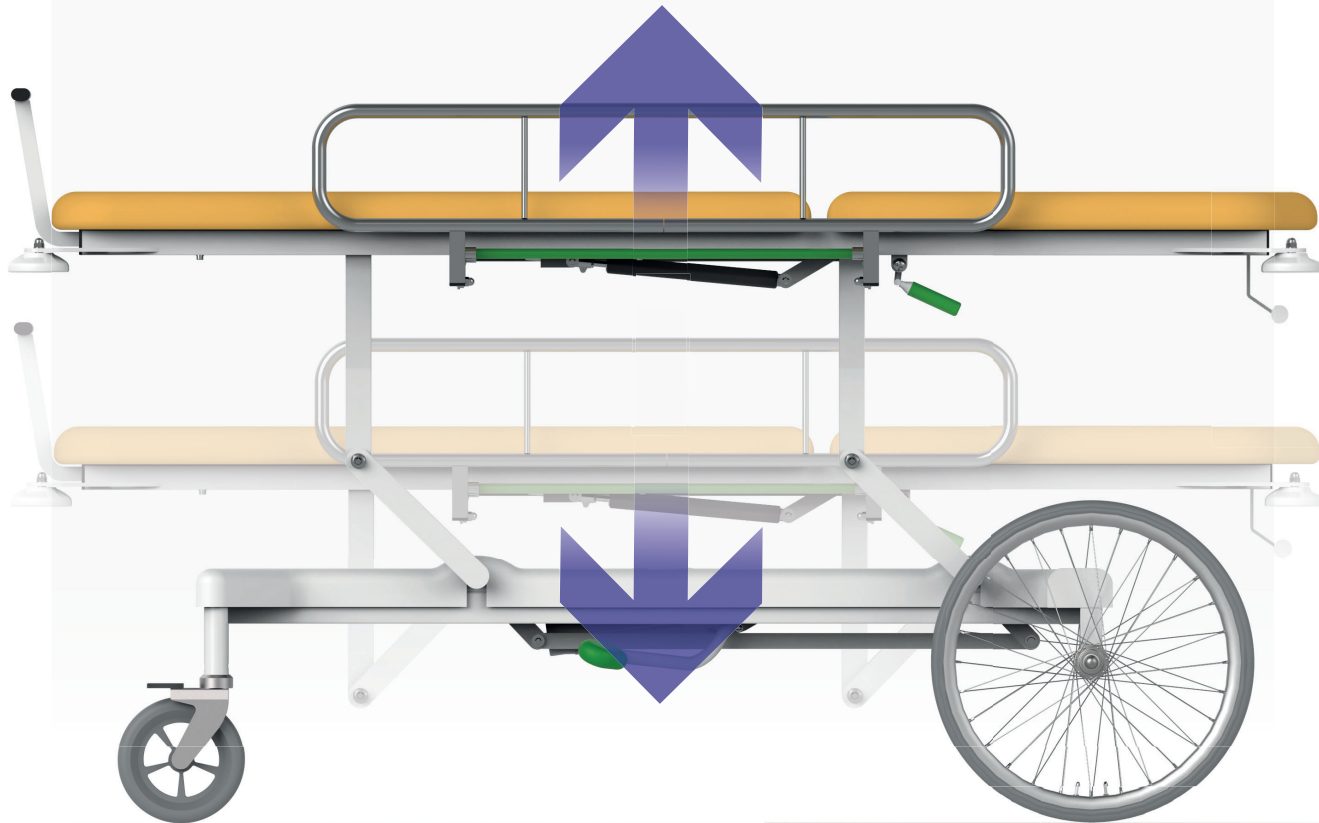



Comfortable positioning

Back positioning is aided by a gas piston rod that allows instant adjustment to the desired position without physical effort.

Safe side rails


The protective side rails along the mattress platform can be easily folded down to provide immediate access to the patient.

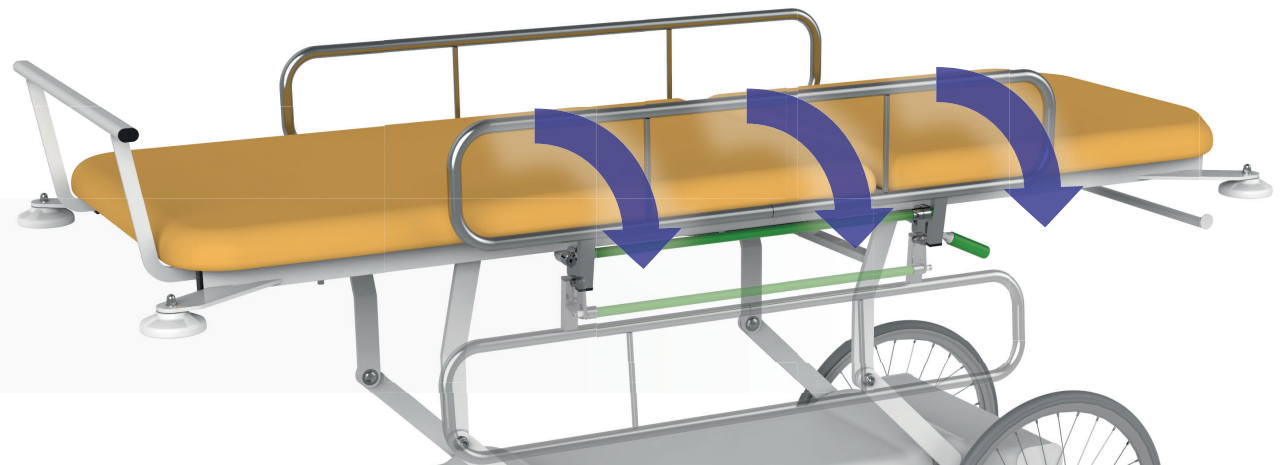

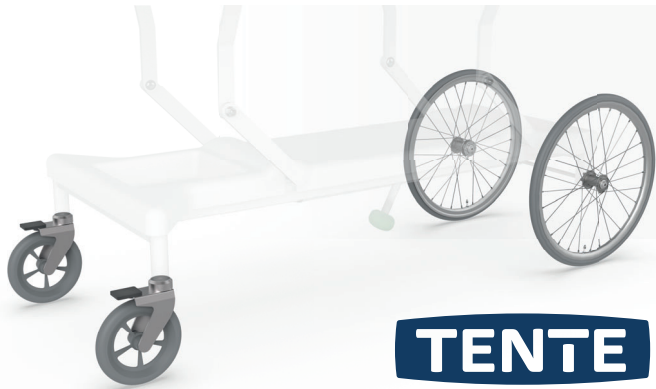


Castors

TENTE's high quality medical castors reduce the physical strain on nursing staff during transport and handling. The castors are made of permanent, non-marking material.

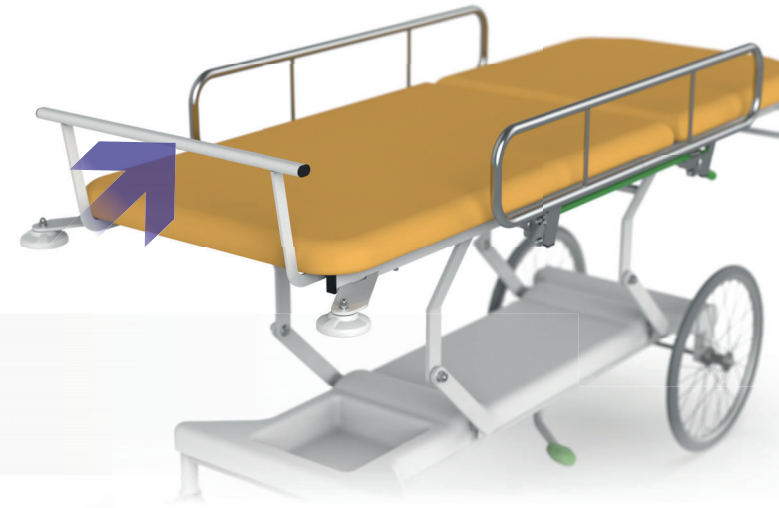


Transport handle in the foot section

For easy patient transport, the trolley is equipped with a handle in the foot section.

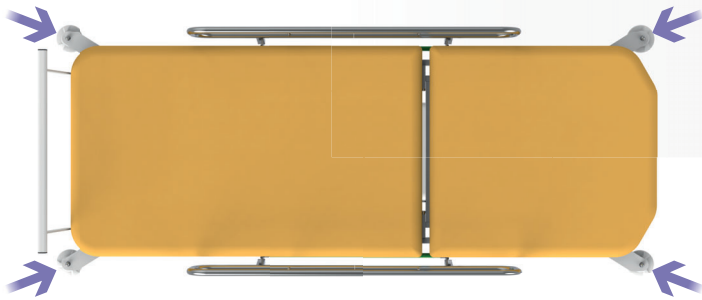

Bumper wheels

The impact castors located in the corners of the mattress platform are a great help especially during transport.

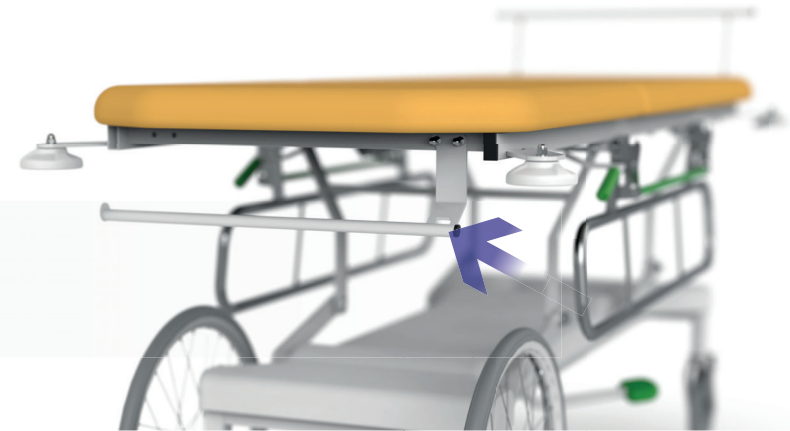

Paper roll holder

Holder underneath the mattress platform for holding hygiene paper.


Storage space


The personal belongings or necessary accessories can be stored in the space under the mattress platform.



Plastic undercarriage cover

The plastic cover makes daily cleaning of the trolley easier and faster. It also protects the chassis mechanisms, thus extending the life of the individual components.

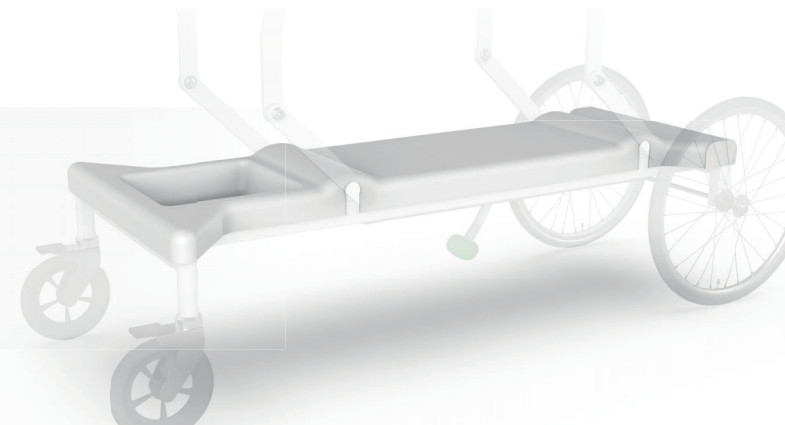

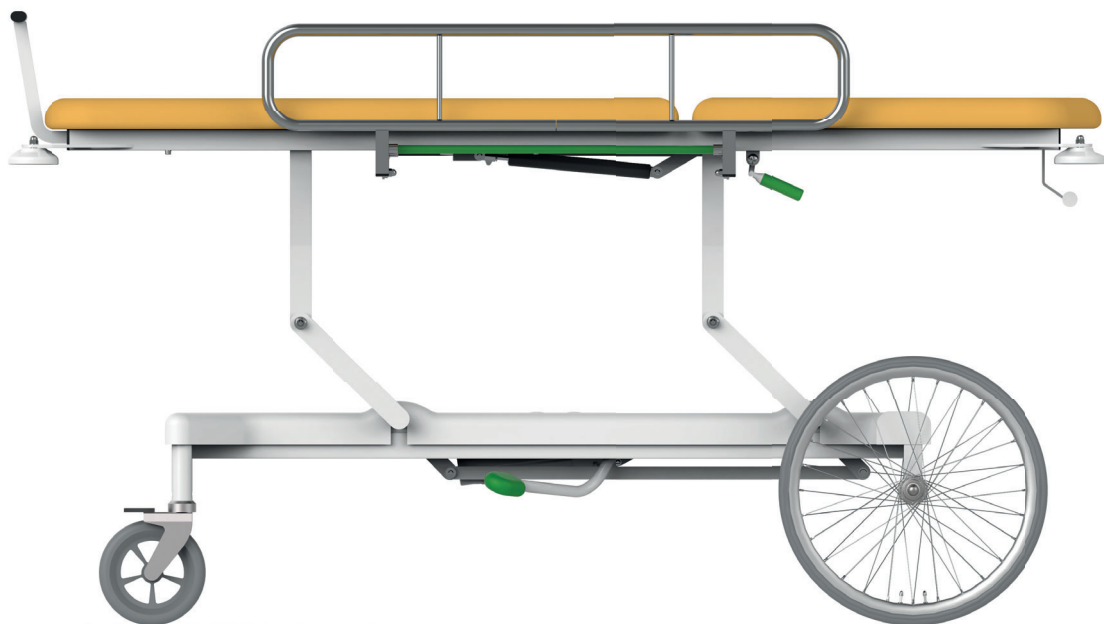

Technical parameters

Outer bed dimensions	208 x 88 cm	Type of siderails.....	Tiltable, stainless steel
Mattress platform	200 x 72 cm	Castors	Metal, bantam, wire
Mattress platform height.....	57 - 94 cm	Material of mattress platform.....	Medical leatherette
Back part positioning	0 - 85°	Plastic undercarriage cover.....	Yes (optional)
Safe operating load.....	230 kg	Paper roll holder.....	Yes (optional)
Maximum patient weight.....	195 kg	Transport handrail	Yes (optional)
Weight.....	90 kg	Bumper wheels	Yes (optional)

