





GRANO

A bed for superior nursing care

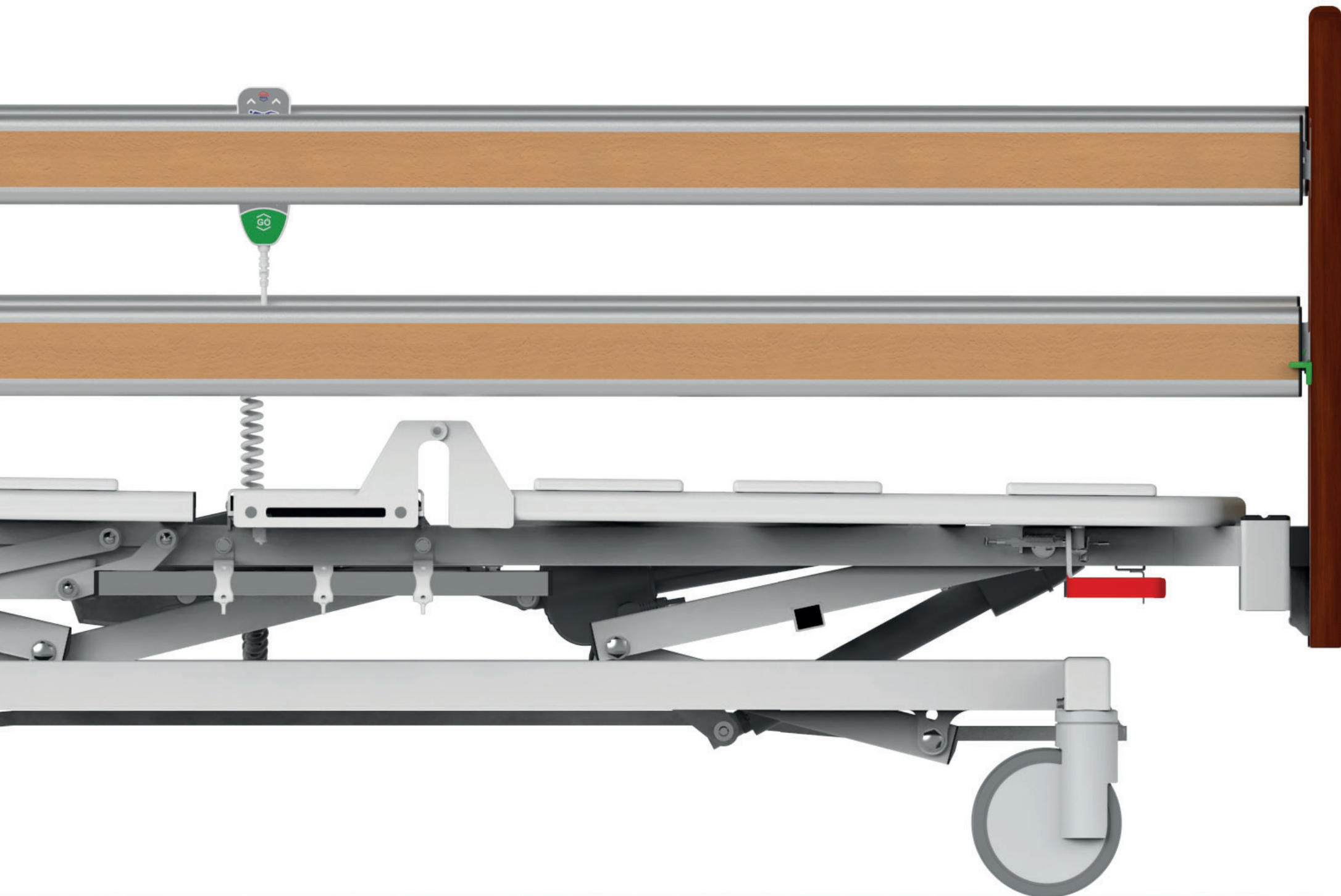

GRANO


GRANO beds are designed for clients of care facilities, such as homes for the elderly, hospices, retirement homes and hospital wards with higher demands on comfort and functionality. The automatic brake system, innovative positioning mechanisms and the patient controller with integrated lamp fully meet the high requirements for safety and comfort of clients. The beds are equipped with electric positioning of the back and foot part of the mattress platform.

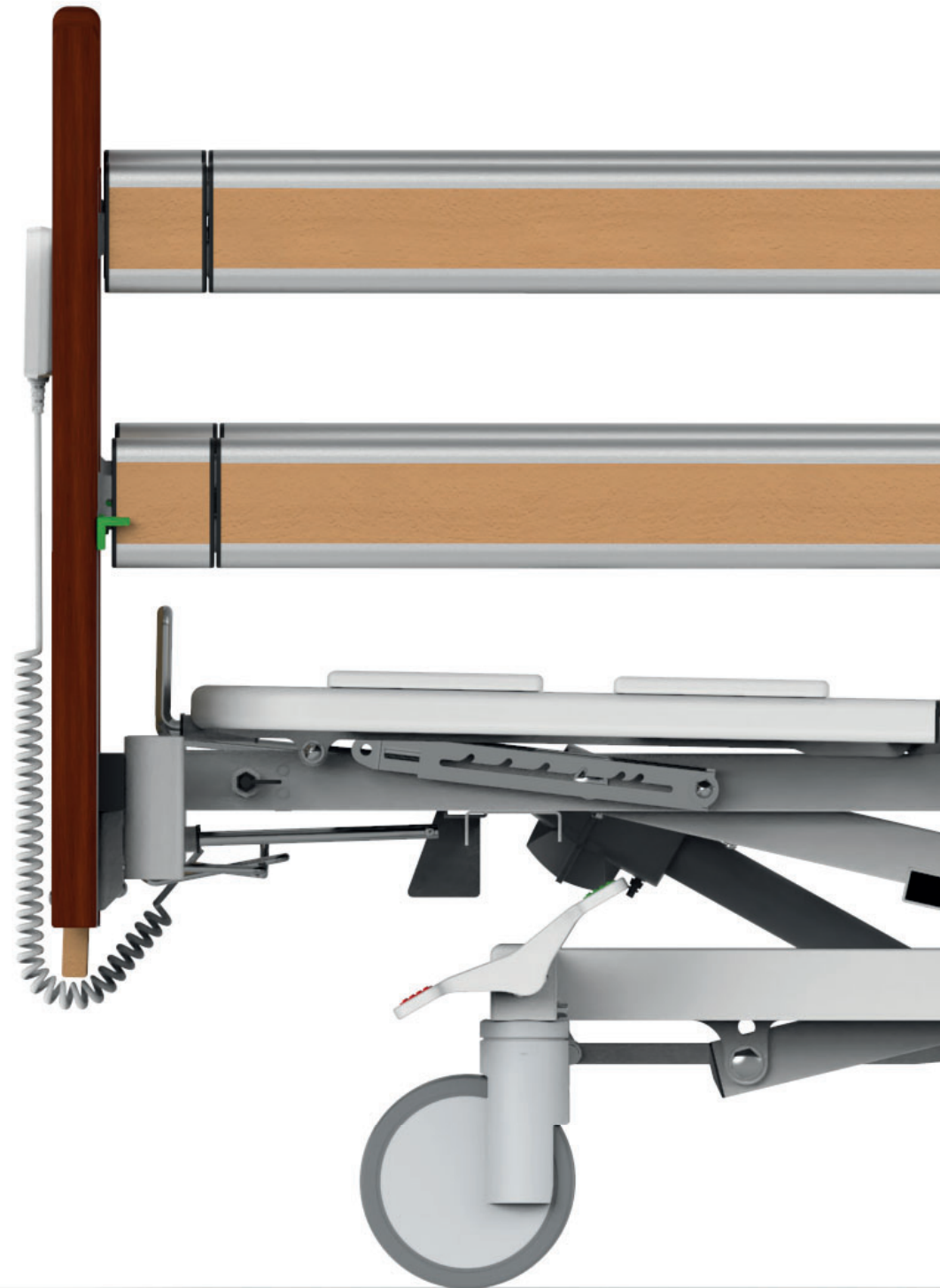




Increased load capacity

The solid and stable construction with a safe operating load of up to 270 kg provides safe support even when caring for patients with higher weight. Thanks to the materials and technologies used, extremely high durability is ensured even under repeated increased loads.

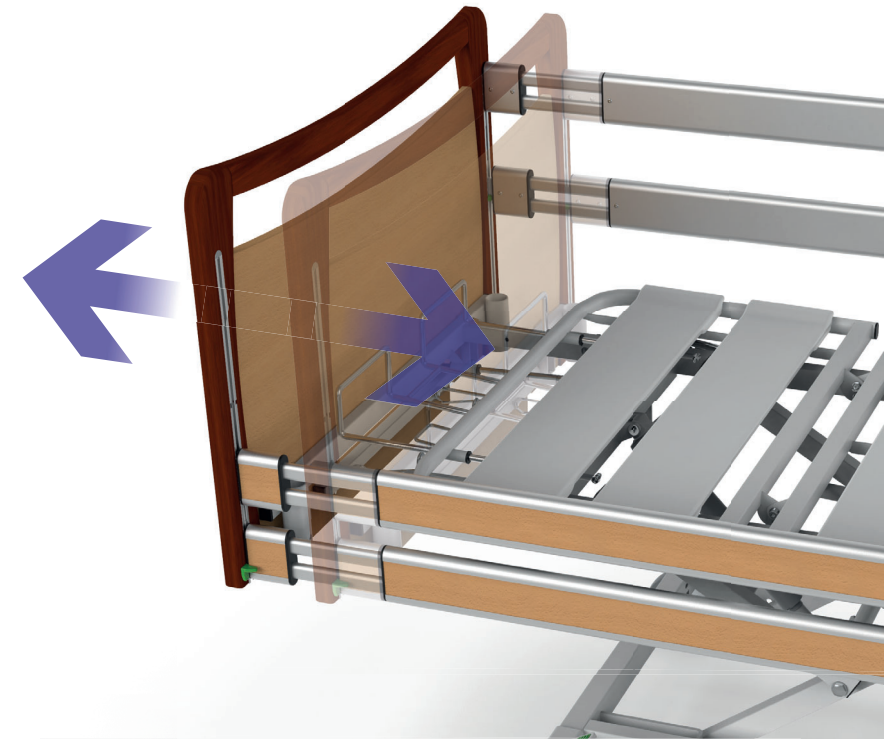


Extension of the mattress platform

Extending the mattress platform increases comfort even for taller patients and makes it easier for staff to handle the patient.

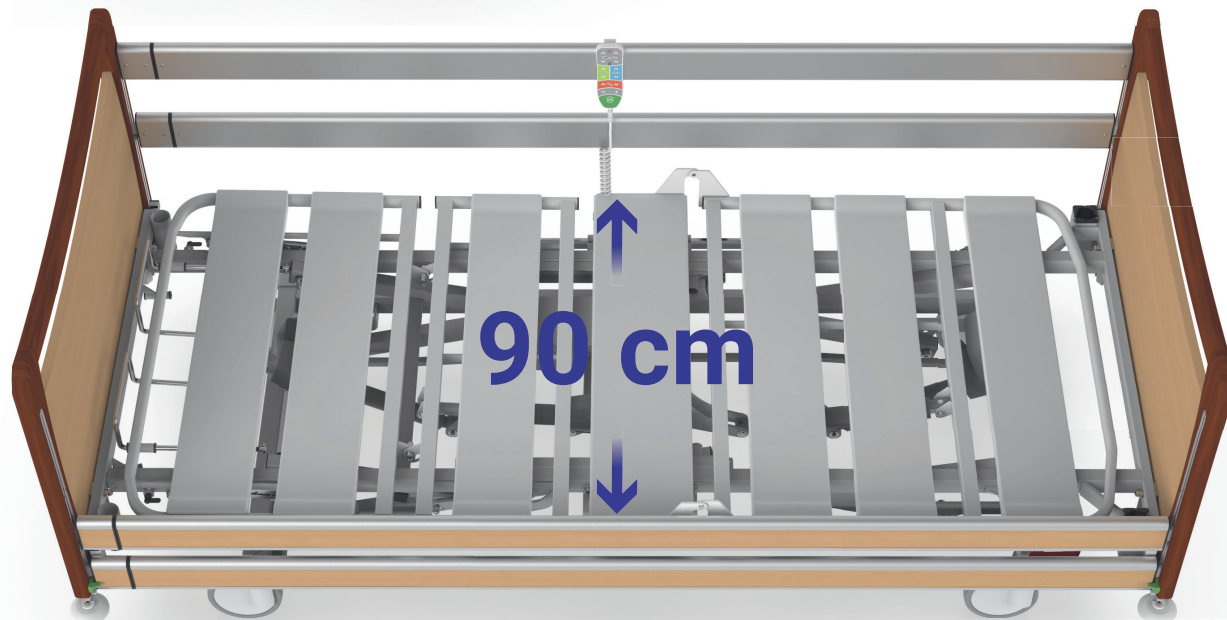


Width of the mattress platform

The optimal width of the mattress platform ensures maximum comfort for the patient and easy handling of the bed.

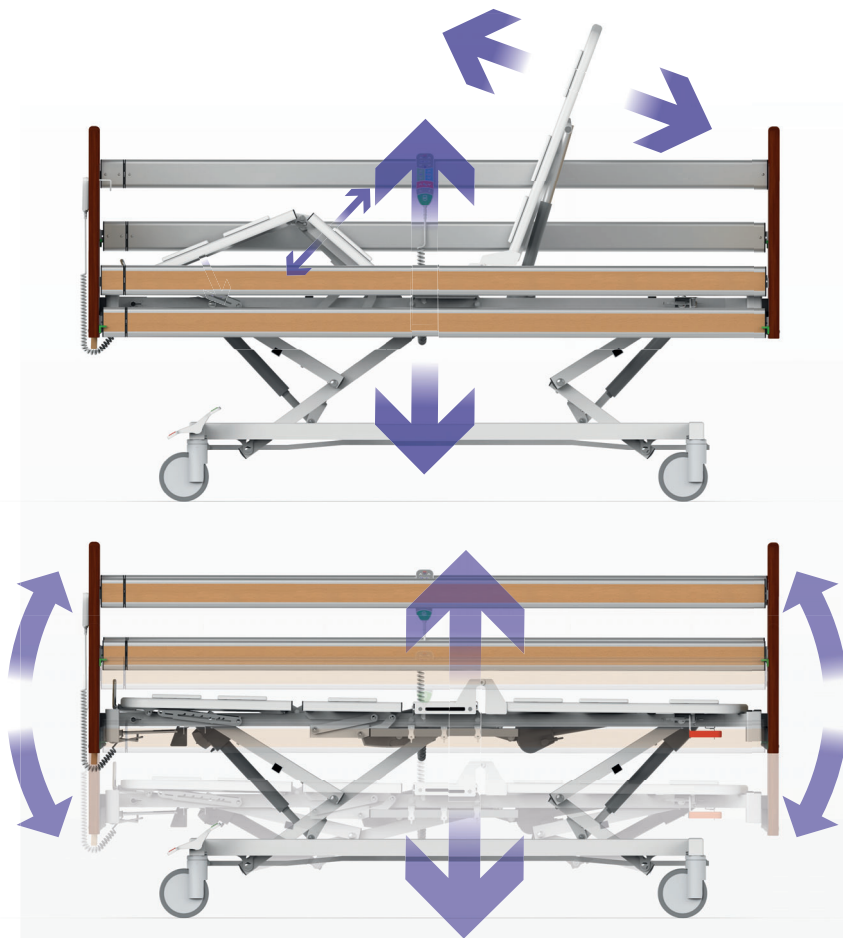

Comfortable positioning


The increased range of the bed lift reduces the strain on staff and facilitates the client handling. The reduced minimum height of the mattress platform assists the client in exiting the bed. The electric positioning of the back and thigh section in combination with the mechanical positioning of the calf section allows a wide range of adjustment of the individual parts of the mattress platform. The staff and client can conveniently adjust the individual positions using the hand control.

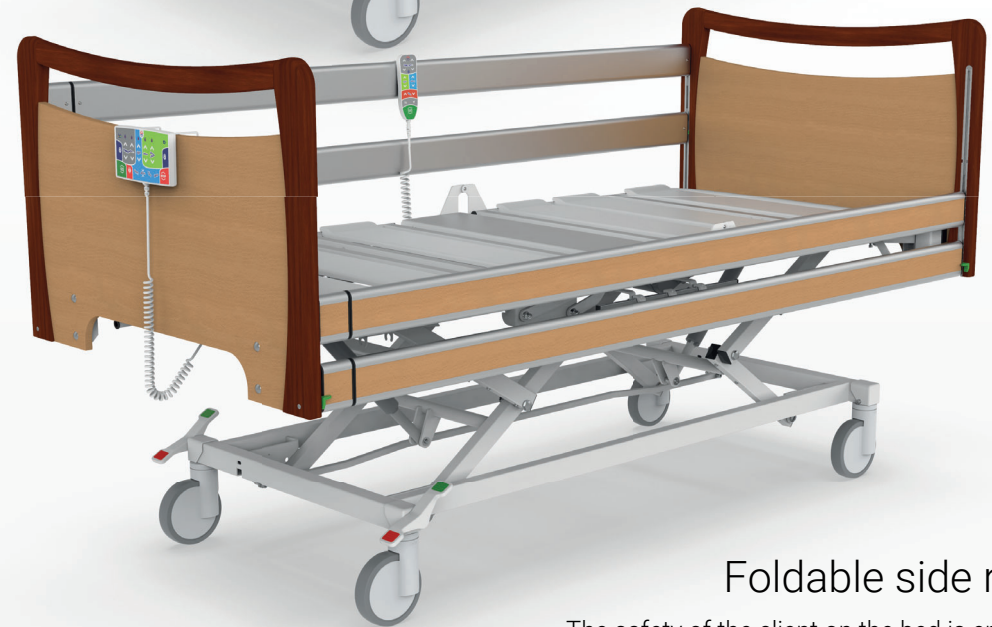


Foldable side rails

The safety of the client on the bed is ensured by folding side rails. The easy folding of the side rails with two independent movements significantly reduces the risk of accidental lowering. The folded side rails are positioned sufficiently below the mattress level to allow easy exit from the bed.



CPR

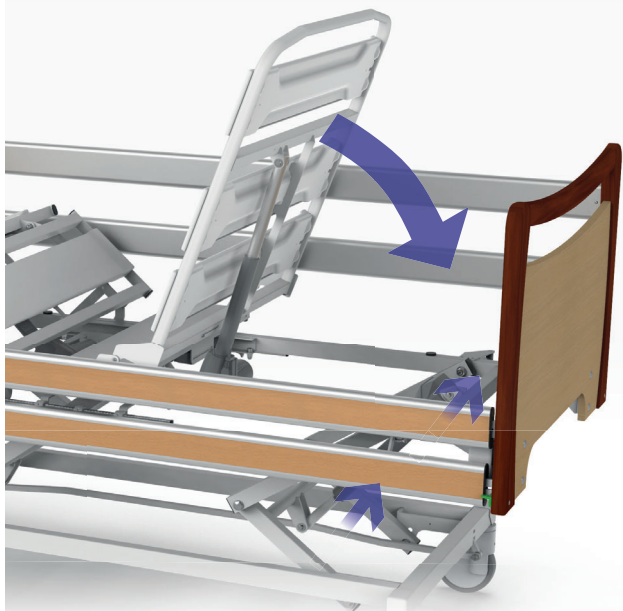

Universal accessories holders

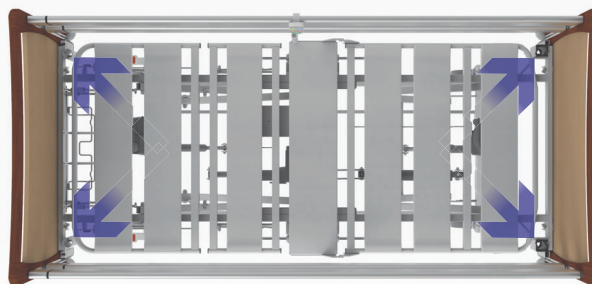


Bumper wheels

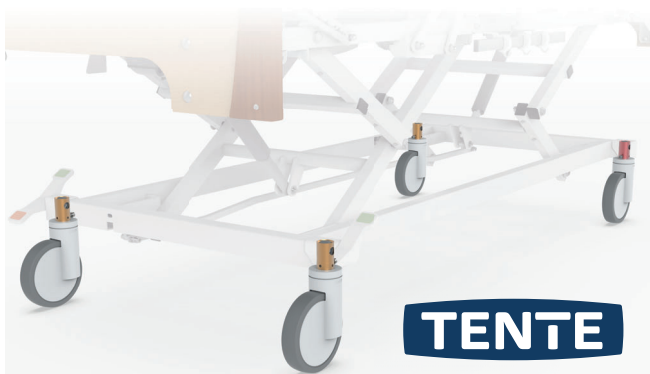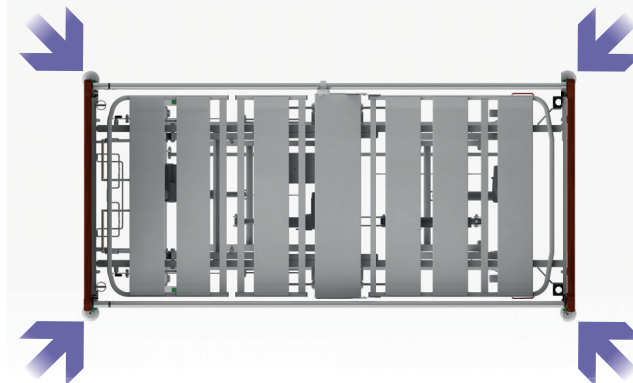

Plastic undercarriage cover





Plastic wheels with central brake

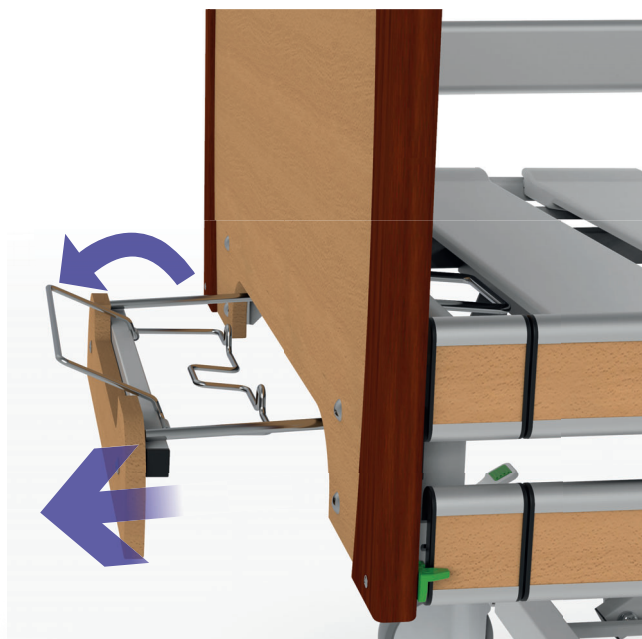

Bed linen shelf

Removable plastic lamellae

The high-quality, easily washable mattress platform made of plastic removable slats ensures high strength and durability even when roughly handled.

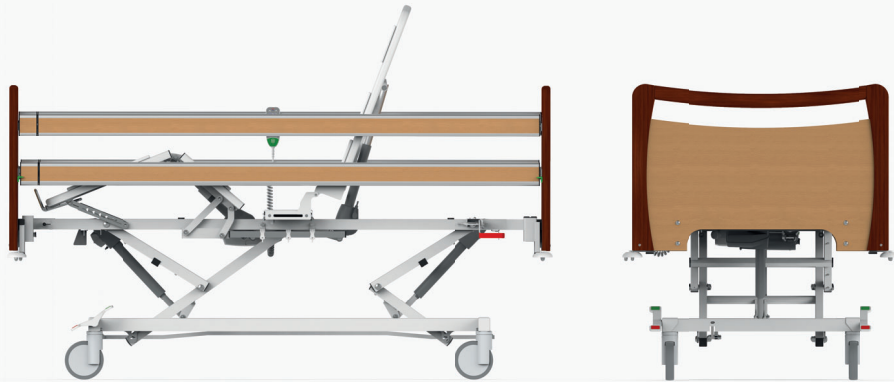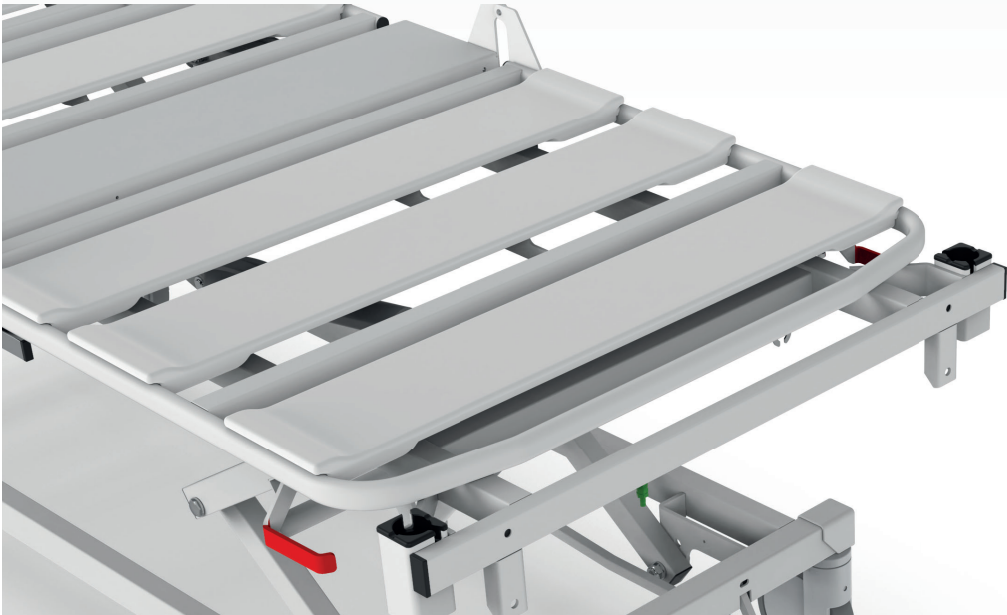


Technical parameters

Mattress platform (length x width)	200 x 90 cm
Outer bed dimensions (length x width)	210 - 213 x 102 - 103 cm
Mattress platform height	40 - 85 cm (electrical)
Mattress platform extension	20 cm (optional)
Calf part extension	Yes
Weight of the bed	125 kg
Safe operating load of the bed	270 kg
Maximum patient weight	235 kg
Application environment for bed use	2, 3
Mattress platform parts - removable plastic lamellae	Yes
Trendelenburg / Antitrendelenburg	16° / 16° (electrical)

Choice of design and decor

Being able to choose the design and colour decor of care beds allows you to create an esthetic and welcoming space that can help improve mood, comfort and reduce stress.

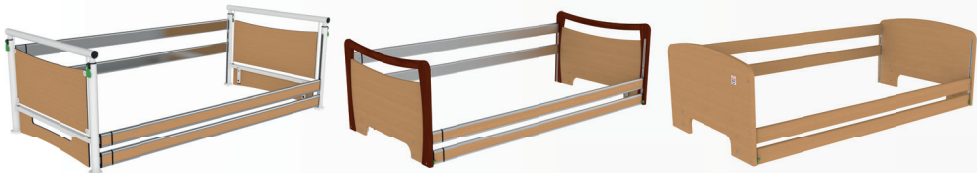

Mattress

Recommended mattress (length x width)	200 x 90 cm
Mattress thickness in normal using	14 cm
Recommended minimal mattress size	12 cm
Recommended maximal mattress size	22 cm

