

LAGUNA

Shower bed for decent and reliable daily hygiene for immobile patients

Hospital staff mainly appreciate the comfortable and safe personal hygiene care of immobile patients or patients with reduced mobility. The shower couch is particularly suitable for use in hospitals, residential care facilities and homes for the elderly.



Increased load capacity

The solid and stable construction with a safe operating load of up to 230 kg provides sufficient support even when caring for patients with higher body weight.






Side rails

Stainless steel hinged side rails with a special antibacterial coating provide maximum patient safety. Thanks to the large folding angle of the side rails, they allow good access to patients.

Comfortable positioning

The shower bed is hydraulically adjustable and completely independent from power source.

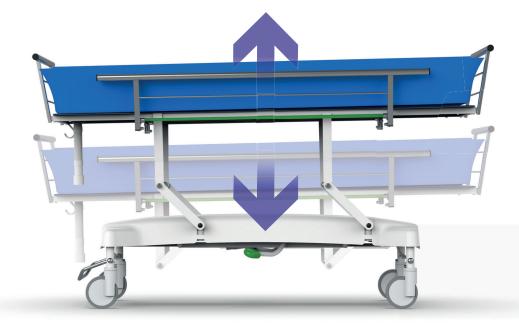



TENTE branded castors

TENTE's high quality anti-static medical castors reduce the physical strain on nursing staff during transport and handling. The castors are made of permanent non-marking material. Castor covers facilitate cleaning, protect bearings from dirt and significantly extend their durability.

> Central brake A central brake operated by a metal bar ensures safe handling of the shower bed.



Because safety, dignity and comfort are most important



Shower tub

The removable shower tub is made of non-slip waterproof material. The bedding surface of the shower tub consists of a one-piece surface with longitudinal drainage grooves. A removable cushion is included in the tub for added patient comfort.

Bed surface

The shower bed has a one-piece loading surface made of HPL. The bedding surface is slightly inclined towards the outlet. This makes it easier to drain the remaining water into the flexible drain hose in the foot section of the bed.

Plastic chassis cover

The plastic cover protects the mechanism from the effects of moisture

and also speeds up the daily cleaning of the shower bed.





Technical parameters

Mattress platform (length x width)	
Outer bed dimensions (length x width)	211 x 79,5 cm
Mattress platform height	54 - 91 cm
Safe operating load of the bed	230 kg
Maximum patient weight	195 kg
Weight of the bed	80 kg
Lifting of the mattress platform	Hydraulic
Mobile	Yes
Type of brakes	Central
Length of drain hose	30 - 70 cm
Castors	150 mm, plastic

PROMA REHA®



www.promareha.com

Riegrova 342, 552 03 Česká Skalice, Czech Republic Tel: +420 775 776 122, info@promareha.com

02/2023-a-EN MD