


PROMA REHA®

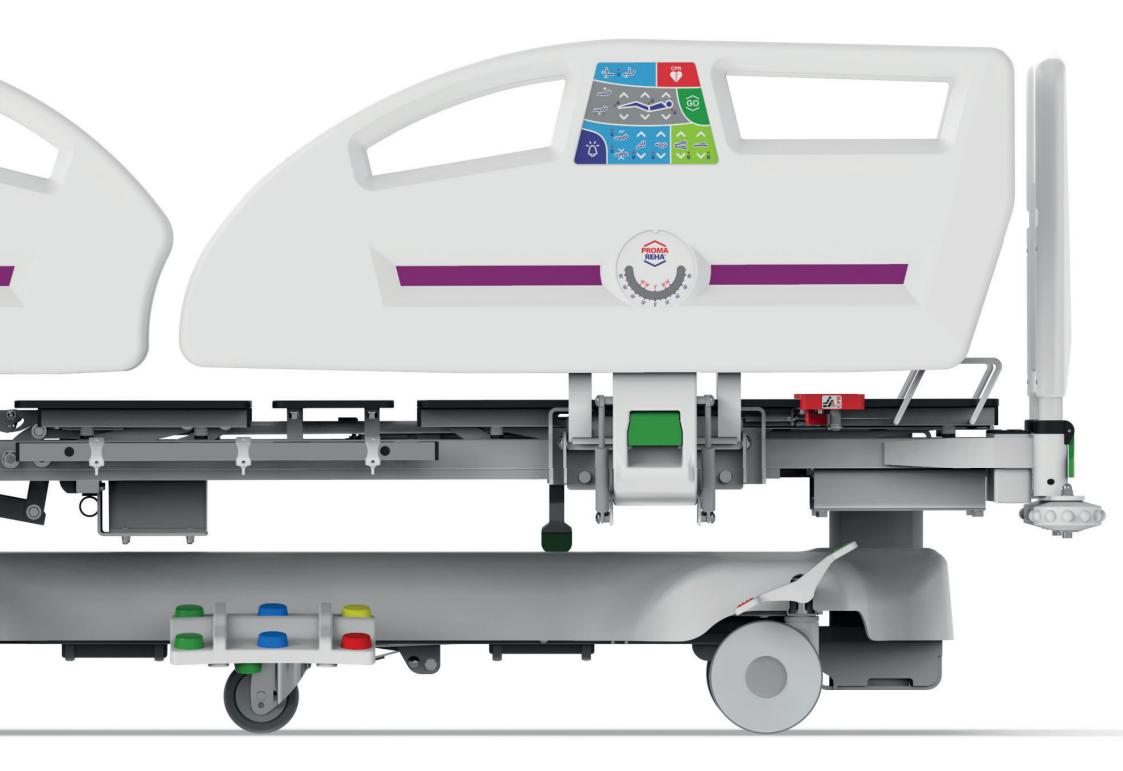


## VISTA

# The most advanced bed for critical and intensive care

The fully electrically adjustable VISTA bed is a huge help for medical staff. Thanks to modern technologies, it reduces the time when handling the patient or during some examinations. It offers lateral therapy and bedside x-rays. Hospital staff and patients alike will appreciate the well-thought-out design, easy handling, comfortable and intuitive control and a whole range of preprogrammed bed functions.

Unique colour-coded backlit controls without complicated touchscreens make it very easy and intuitive for staff and patients to operate not only in crisis situations but also for people with impaired vision, without glasses or at night.

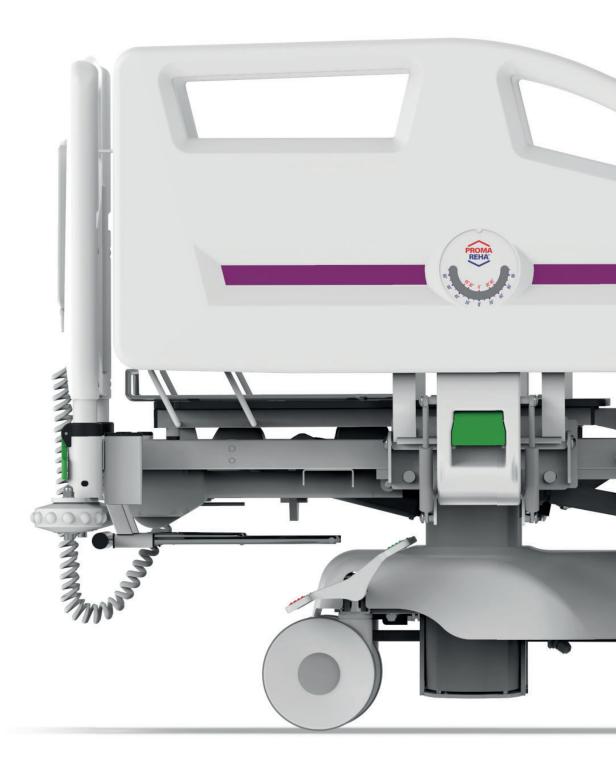


#### Increased load capacity

The solid and stable construction allows a safe load capacity up to 300 kg.



#### X-rayable mattress platform

The design of the mattress platform allows comfortable and safe continuous scanning with C-arm and conventional X-ray equipment.

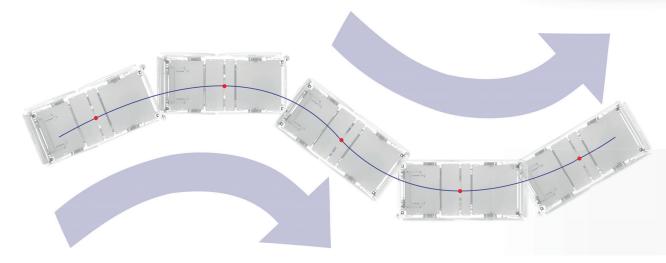

### Fifth wheel

The fifth wheel ensures easy and intuitive handling of the bed. This feature is particularly useful in rugged areas and when the bed is handled by one person.





#### Reversible lateral tilt


Significantly reduces the physical strain on the nursing staff during daily patient care, especially when turning the patient, performing personal hygiene and changing the bed linen.

#### Side rails suitable for higher mattresses

Specially designed siderails with a height of 45 cm above the bed surface guarantee maximum protection against falling out of bed even when using an anti-decubitus mattress up to 23 cm high.

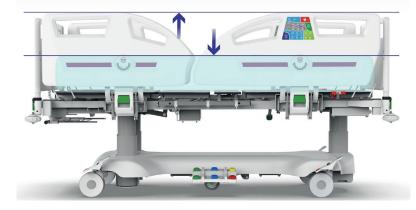

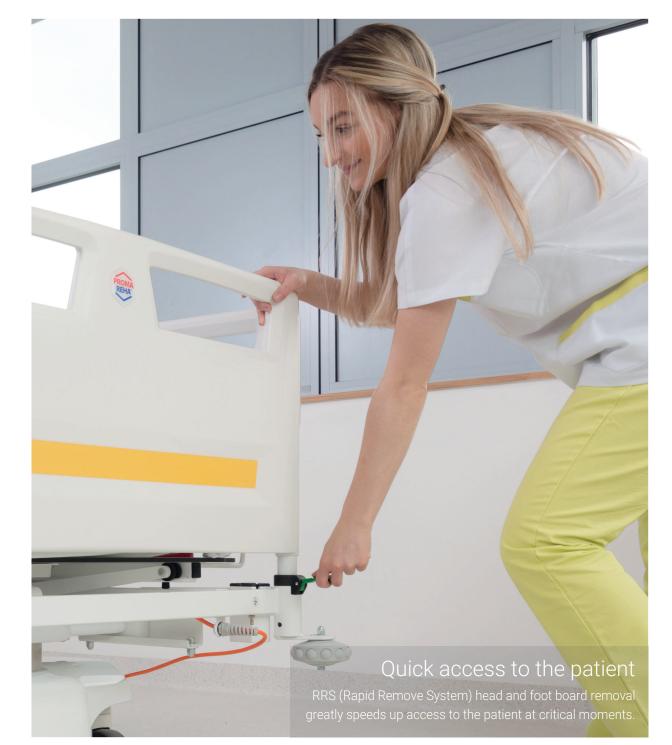

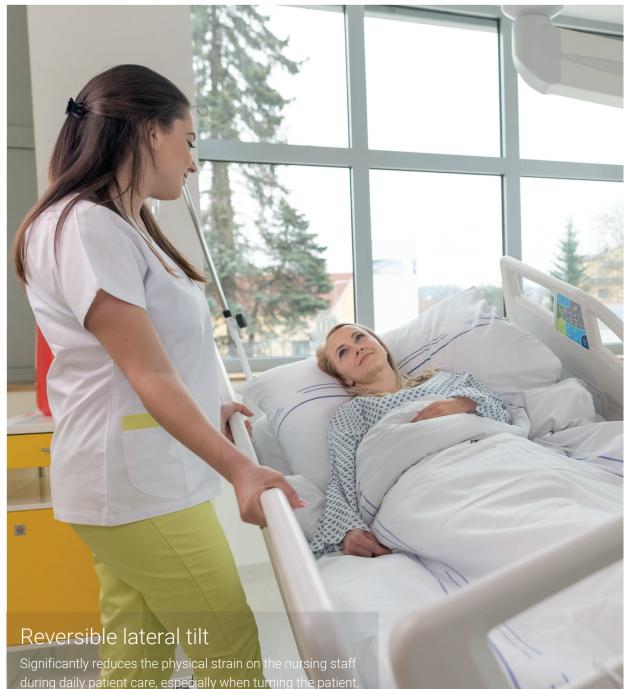

Integrated controls in the side rails

PROMA

With activation button to prevent accidental use

Easily removable boards With safety lock

#### Nurse control with LCD

CPR Electric and mechanical

Backup batteries

Bed linen shelf

Full length plastic side rails For maximum patient safety

during daily patient care, especially when turning the patient, performing personal hygiene and changing the bed linen.



#### Horizontal accessories holder

For hanging additional equipment

#### X-rayable mattress platform

Using both C-arm and conventional X-ray cassettes

#### TENTE brand wheels

150 mm plastic with central brake

Increased load capacity 300 kg


#### Chassis illumination

#### Bed and back tilt indicator

For easy adjustment of the treatment position


## Universal accessories holders

Intuitive controls with colour-coded functions

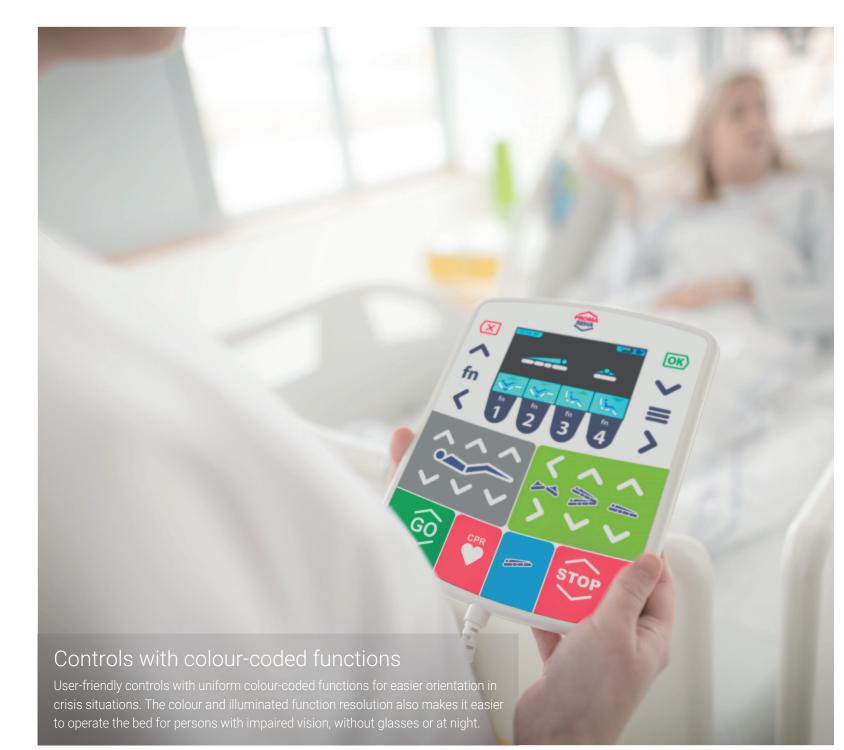





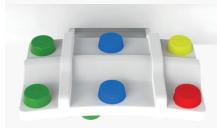





Fowler's position



Height of mattress platform 44 cm - 84 cm



Position of the cardiac chair



Trendelenburg / Antitrendelenburg



An electric bed extension of up to 20 cm brings comfort even to taller patients.

#### Optional equipment



Weighing system Weight is one of the important parameters of patients' health status, especially in acute care.


> Electrically driven fifth wheel Handling the bed has never been so easy as it is now with the electric drive.

#### Automatic lateral therapy

It significantly reduces the risk of postoperative respiratory complications.

### PUSH & GO



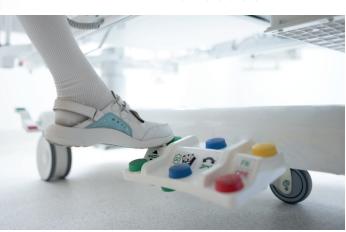





Foot controls

controls. All basic functions are controlled with your foot and your hands are free.

#### Patient control

The bed can be optionally equipped with an ergonomic patient control, which is always at hand of the patient and allows easy control of all basic functions.

> X-ray holder for side mounting For more comfortable handling of the X-ray cassette without the need to position the patient.

Give yourself a free hand with the foot



#### Technical parameters

Mattress platform length Mattress platform width Outer bed dimensions width Outer bed dimensions length Mattress platform height Mattress platform extension Height of headboards and sideboards	90 cm 
above the loading surface	
Back length	
Pelvic length Thigh length	
Calf length	
Weight of the bed	153 kg
Safe operating load of the bed	
Maximum patient weight	235 kg
Mattress thickness in normal using	
Recommended minimal mattress size	
Recommended maximal mattress size	
Recommended mattress (lenght x width)	200 x 90 cm
TR / ATR angle	16° / 14°
Lateral tilt	
Castors	150 mm
	150 mm
Castors Autoregression back part	
Autoregression	14 cm
Autoregression back part thigh part	14 cm
Autoregression back part	14 cm 7,5 cm

Weight range	) - 300 kg
Accuracy of weighing in absolute mode	± 0,5 kg
Accuracy of weighing in differential mode	± 0,5 kg



PROMA REHA®



www.promareha.com

Riegrova 342, 552 03 Česká Skalice, Czech Republic Tel: +420 775 776 122, info@promareha.com

02/2023-a-EN