



# decu® plus combi

width 80/90 cm

load capacity 200 kg




III. deg. pressure ulcers

air

Dynamic anti-decubitus air mattress



The traditional air mattress conception is widely used all over the world in home, social and hospital care. Replacement mattress with height 13 cm air + 5 cm foam for decubitus stage II. to III./medium to very high with load capacity 200 kg. The mattress is composed of 20 cells, 4 cells under the head have a pillow function and 6 cells have a micro-ventilation system. Available in two widths 80 or 90 cm. The safety Quick Release connectors for achieves zero pressure in high risk areas (option). Semi-automatic compressor with automatic pressure settings according to the patient weight.

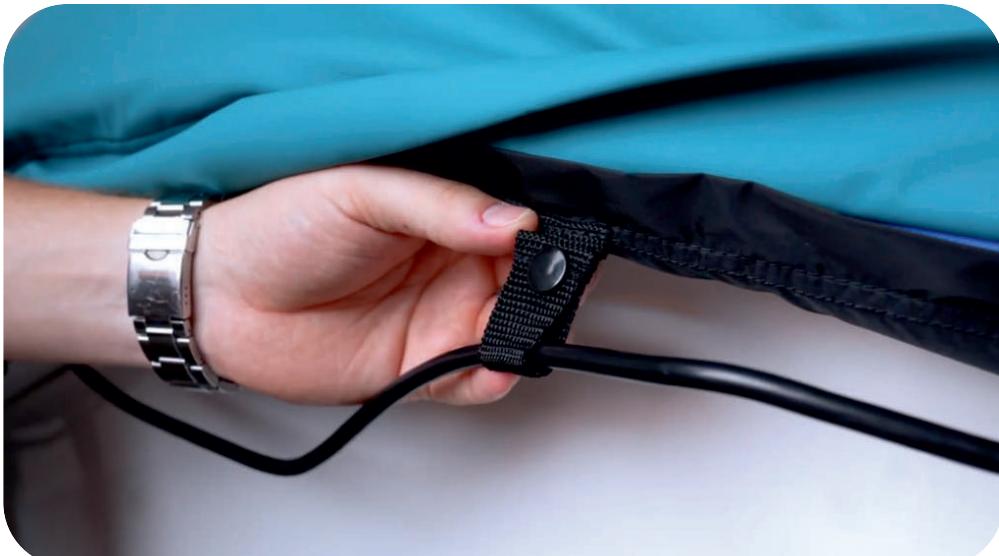

# MATTRESS

- 6 cells under hull with micro-ventilation function to reduce heat and moisture
- 2-way stretch, waterproof, breathable, fire-retardant cover
- Easy maintenance - washable up to 95°C
- Zipper 360° with perimeter protect flap which reduces the risk of contamination
- Integrated power cable fixation system
- 4 straps for easy fixation of the mattress to the bed
- Quick release CPR valve for acute resuscitation
- QR - safety Quick Release connectors for easy disconnection of the cells and to achieve zero pressure in high risk areas to reduce the risk of decubitus - 4 / 20 pes (option)

## TECHNICAL PARAMETERS

<b>name</b>	decu plus combi
<b>type</b>	air mattress
<b>stage of decubitus</b>	II. - III. (medium to high)
<b>generic group</b>	dynamic anti-decubitus air mattress
<b>Basis UDI-DI</b>	8594218991DAAME3
<b>REF / GTIN (90 cm)</b>	1-329-01 / 08594218993299
<b>REF / GTIN (80 cm)</b>	1-328-01 / 08594218993282
<b>package size</b>	48 x 38 x 26 cm
<b>package weight</b>	10 kg
<b>package includes</b>	mattress with cover, compressor, user manual, transport bag
<b>warranty</b>	2 years

# COMPRESSOR

- Automatic pressure setting according to the selected patient weight
- Alternate mode: APT - Alternate Pressure Therapy mode A-B-C (5/10/20 min)
- Static mode: CLP - Constant Low Pressure mode in all cells
- Auto-firm mode: LMP - time-Limited Maximum Pressure mode (30 minutes)
- Transport mode: TM - Transport Mode with memory function and easy disconnection
- Automatic key lock for unintentional changes to compressor settings
- Acoustic and visual alarm: low pressure, power failure, service
- AVS - Anti Vibration System for vibration reduction and quiet operation (noise < 23,8 dB)
- START/STOP system for energy reduction
- Memory function for keeping of the set values
- Integrated compressor holder

## MATTRESS

<b>type</b>	replacement mattress
<b>load capacity</b>	200 kg
<b>height</b>	13 cm
<b>dimension</b>	80/90 x 200 cm
<b>number of cells</b>	20 pes - width 10 cm
<b>pillow function</b>	4 head cell without alternation
<b>micro-ventilation</b>	6 cells under hull
<b>CPR valve</b>	quick release - manual
<b>next functions</b>	integrated power cable fixation system 4 straps for easy fixation to the bed
<b>material</b>	polyurethane
<b>cleaning</b>	thermal disinfection
<b>weight</b>	5,9 kg
<b>cover</b>	PU, 2-way stretch, waterproof, breathable, fire- power supply retardant, 95°C washable, zipper 360°with flap

## COMPRESSOR

<b>type</b>	semi-automatic with LCD
<b>modes</b>	4 - alternate/static/ auto-firm/transport
<b>alternate mode</b>	3-chamber air exchange system A-B-C
<b>time interval</b>	5/10/20 minutes
<b>alarm</b>	acoustic and visual
<b>alarm type</b>	low pressure/power failure/service
<b>key lock</b>	automatic
<b>next functions</b>	memory, Start/Stop, AVS
<b>noise</b>	<23,8 dB
<b>holder</b>	integrated
<b>flow rate</b>	8 l/min
<b>operating pressure</b>	25-86 mmHg
<b>material</b>	ABS housing
<b>IP</b>	42
<b>size</b>	17x26x12,5 cm
<b>weight</b>	3,4 kg

